





Patient App Set-Up MedM

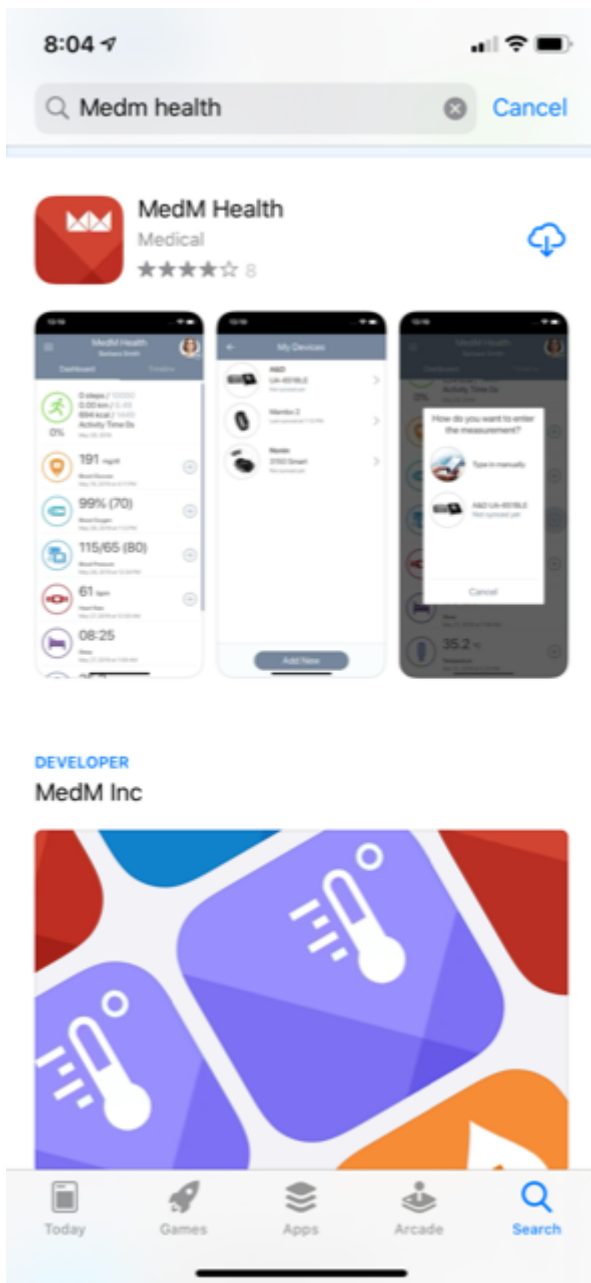

Search and Download MedM Health App

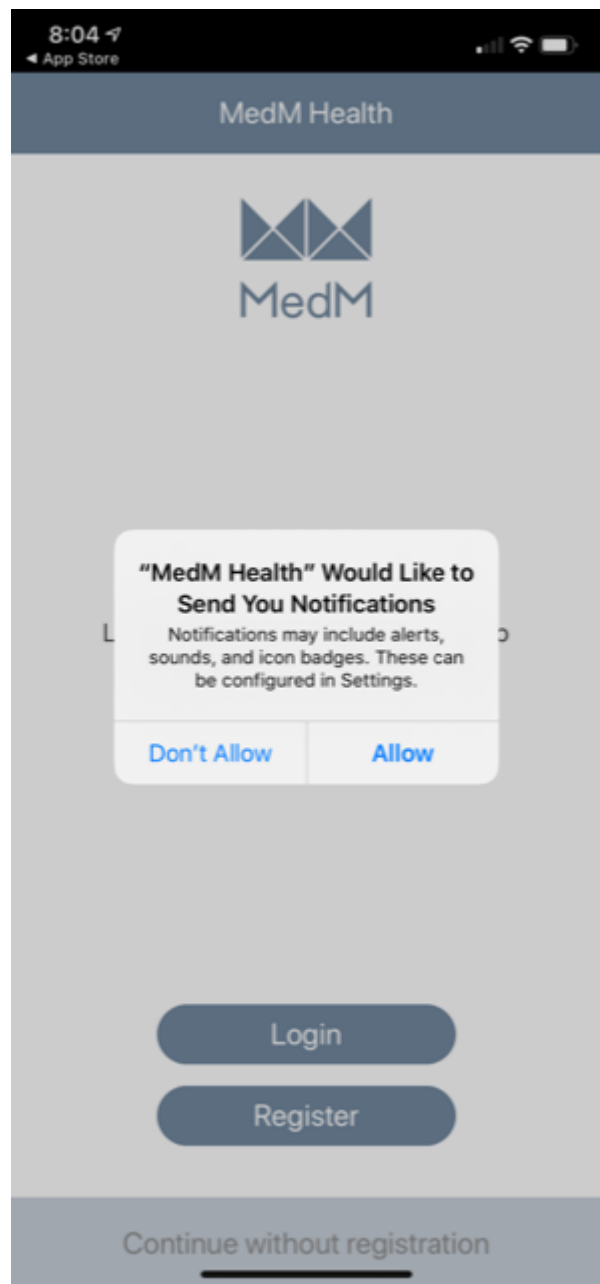

Allow for notifications



Patient App Set-Up MedM




Login or Register with MedM cloud to
get sync and backup your data

Login

Register

Continue without registration

Hit the Register button

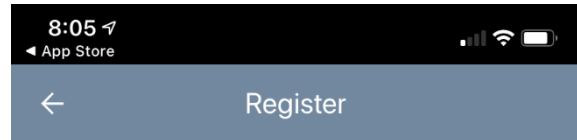

→ robkent97@gmail.com

→ floridarpn



Use 8 or more characters with a mix of letters,
numbers & symbols



I am 16 or older and I agree to the
Terms of Service and Privacy Policy

Register


→ Enter Enterprise Access Code

Enter email and password. Select that
you are 16 or older and agree to terms
and conditions.


DO NOT HIT REGISTER! Select “Enter
Enterprise Access Code” at the bottom.



Patient App Set-Up MedM

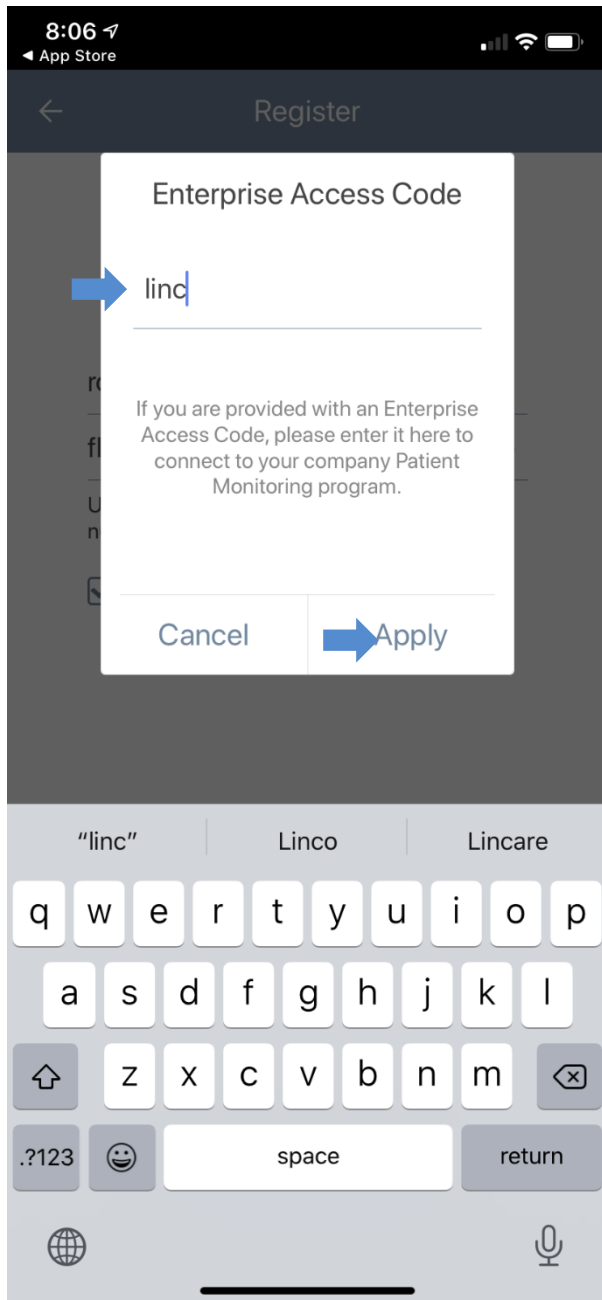

Enter "linc" as the Enterprise Access Code and hit Apply

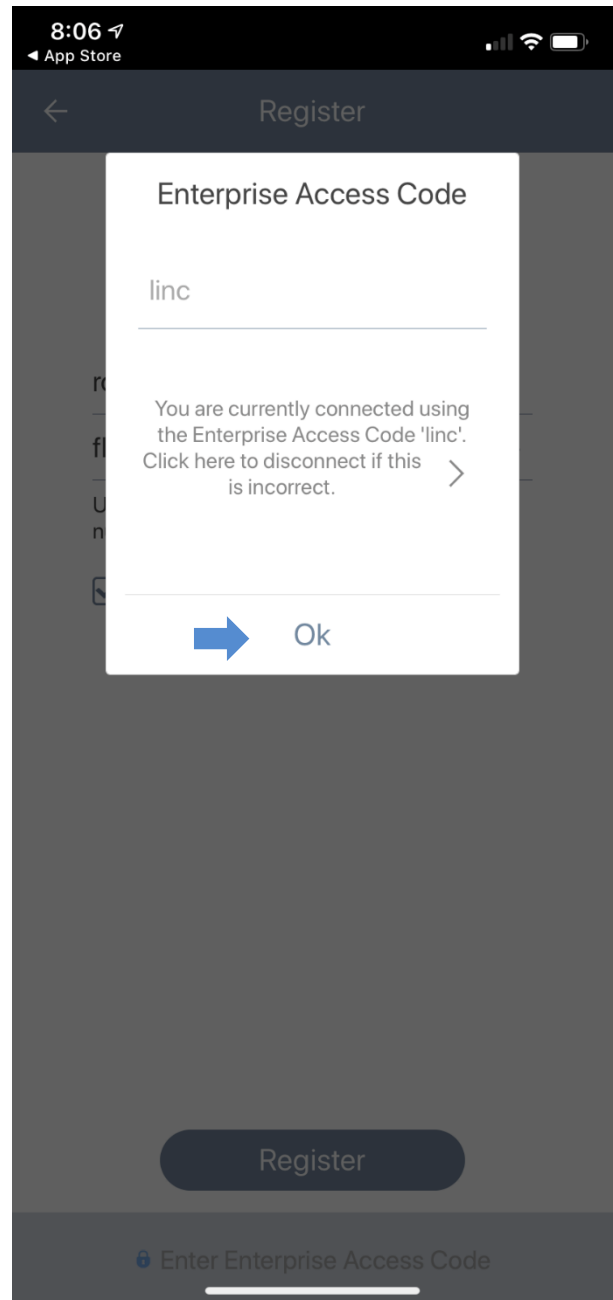


Verify message that you are connected and hit "OK"




Patient App Set-Up MedM

8:06 7
App Store

← Register



MedM



robkent97@gmail.com

floridarpm 

Use 8 or more characters with a mix of letters, numbers & symbols

☒ I am 16 or older and I agree to the Terms of Service and Privacy Policy

 Register


  linc


Verify “linc” in the bottom bar and hit the “Register” button.

8:07 7
App Store

MedM Health
robkent97@gmail.com


Dashboard Timeline

 Blood Glucose
No data yet

 Blood Pressure
No data yet

PIN Protection


Would you like to set PIN protection?

No  Yes


Select “Yes” for PIN Protection and set a 4-digit PIN number. Also set biometric data if prompted and desired.



Open menu at top

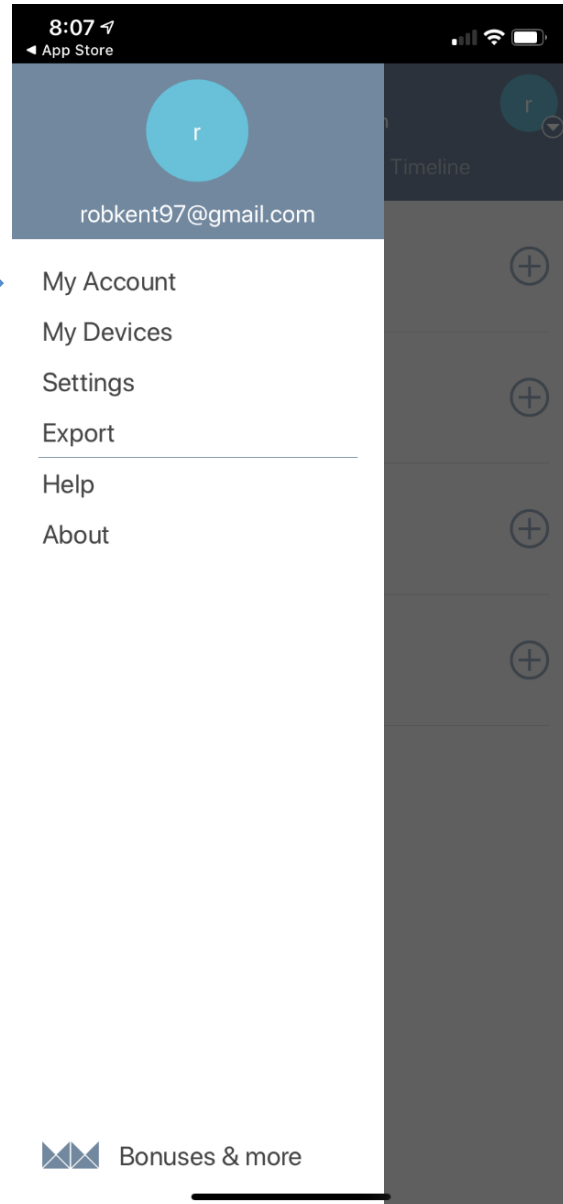

Select My Account

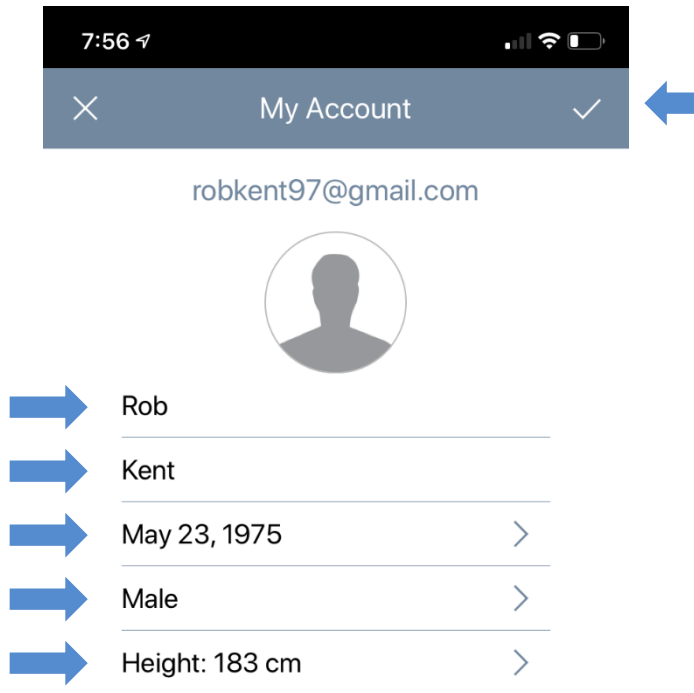

7:56 7:56 7:56

My Account

robkent97@gmail.com

robkent97@gmail.com

Rob

Kent

May 23, 1975

Male

Height: 183 cm

Fill out demographic and height
data and hit check button at top




Jumper Medical Pulse Oximeter




Now attach and turn on your pulse oximeter. It will start streaming readings and will start-up in pair mode.




Open menu at top

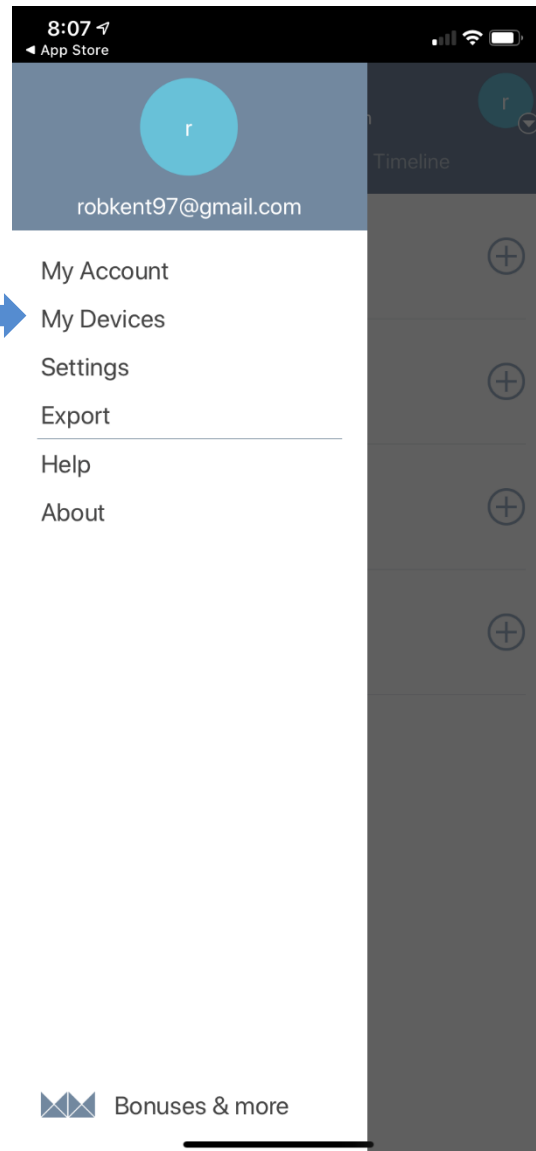

Select My Devices

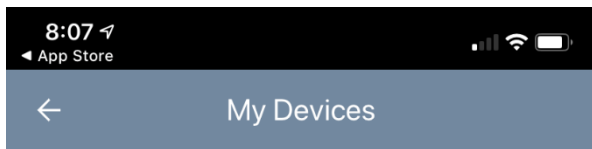

You have no devices. To add a device tap 'Add New' button.

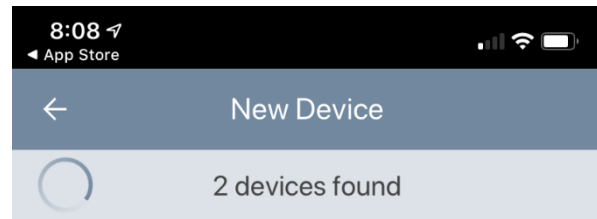

Internal Step Counter
Pedometer based on embedded step counter



Jumper Medical JPD-500D
Oximeter
Signal: excellent (-47 dBm)




Select "Add New"




Select the connected Pulse Oximeter




Patient App Set-Up MedM

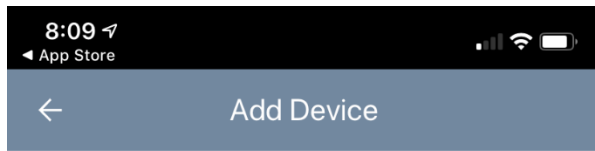

Jumper Medical JPD-500D

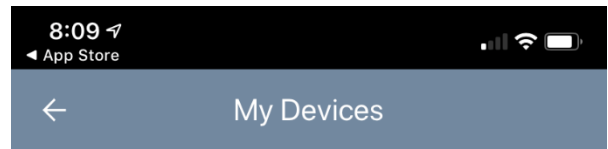


Jumper Medical
JPD-500D
Not synced yet

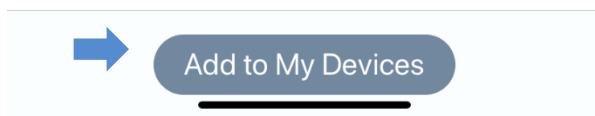


Select “Add to My Devices”




Back out to the Home screen

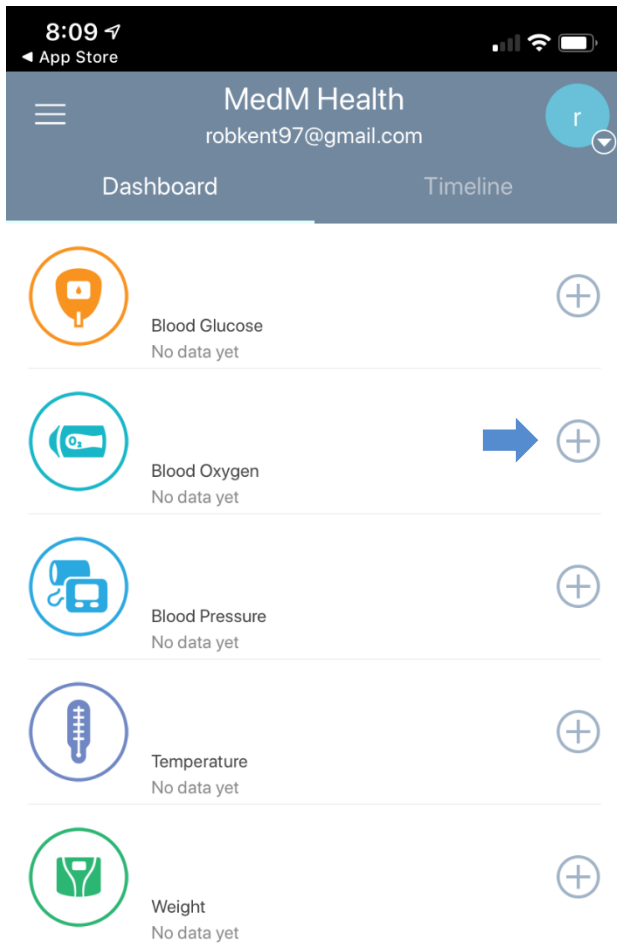

Tap the + sign next to any vital to take a reading

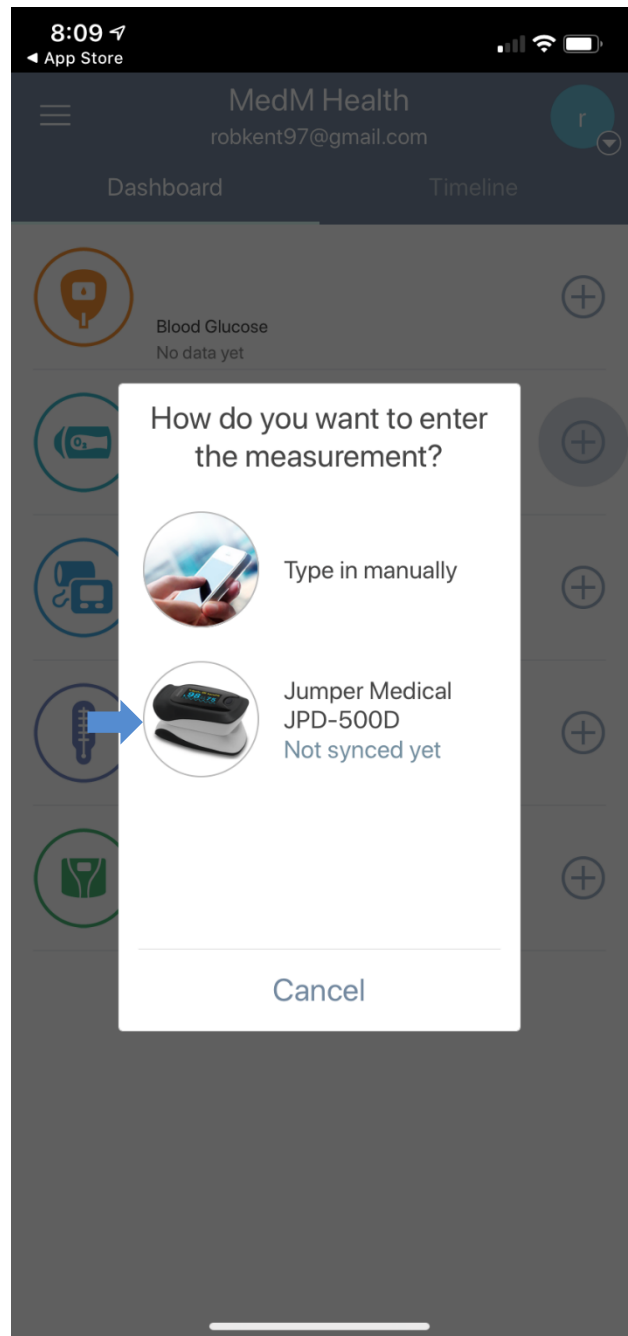

Select the synced device or manual entry



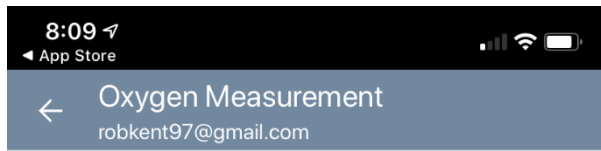


Jumper Medical JPD-500D




Connected. Waiting for data...

Done

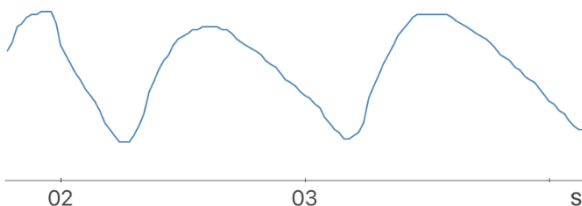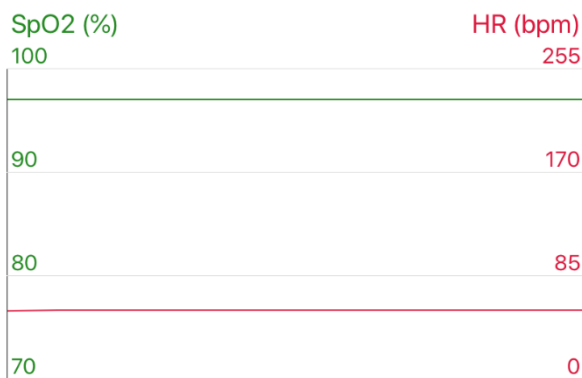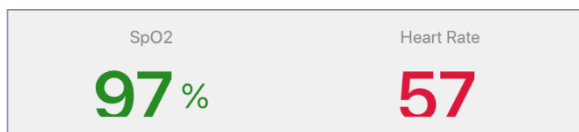

Wait for readings to start streaming. When readings established select "Done"

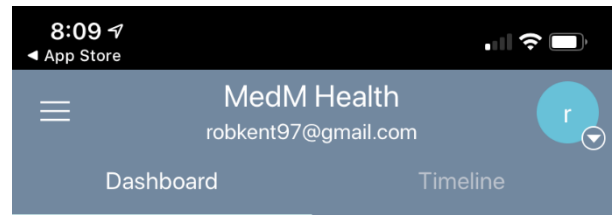

Blood Glucose
No data yet



97% (55)
Blood Oxygen
Sep 28, 2020 at 8:09 AM



Blood Pressure
No data yet



Temperature
No data yet




Weight
No data yet




Latest reading will register on the home-screen and in the Linc database




A&D Blood Pressure Monitor

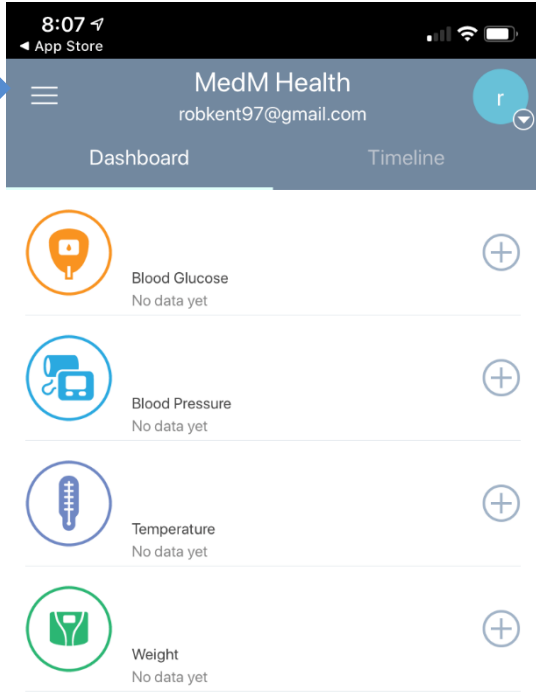


Open menu at top

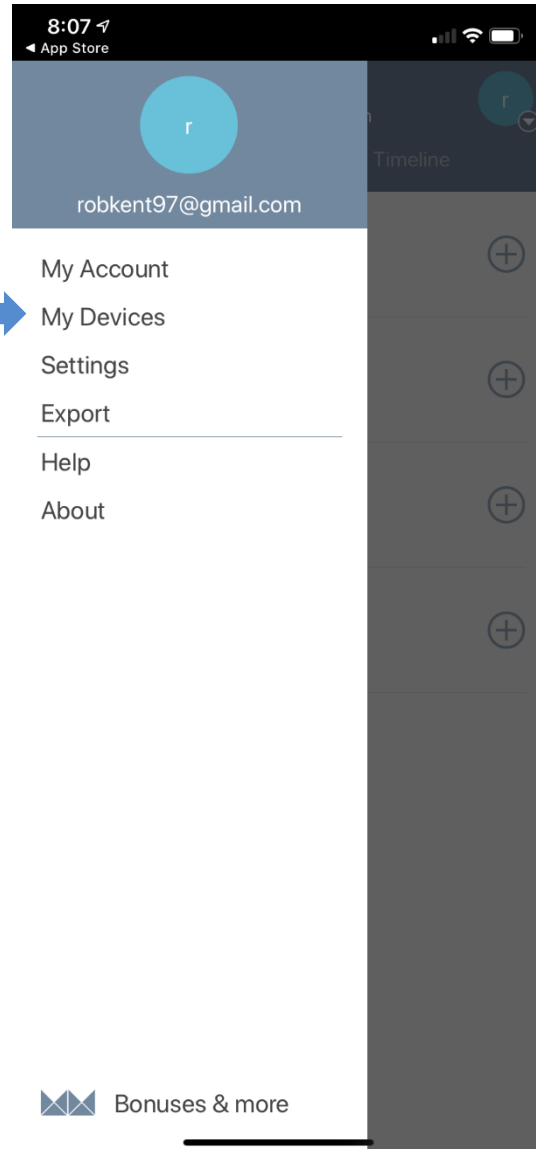


Select My Devices

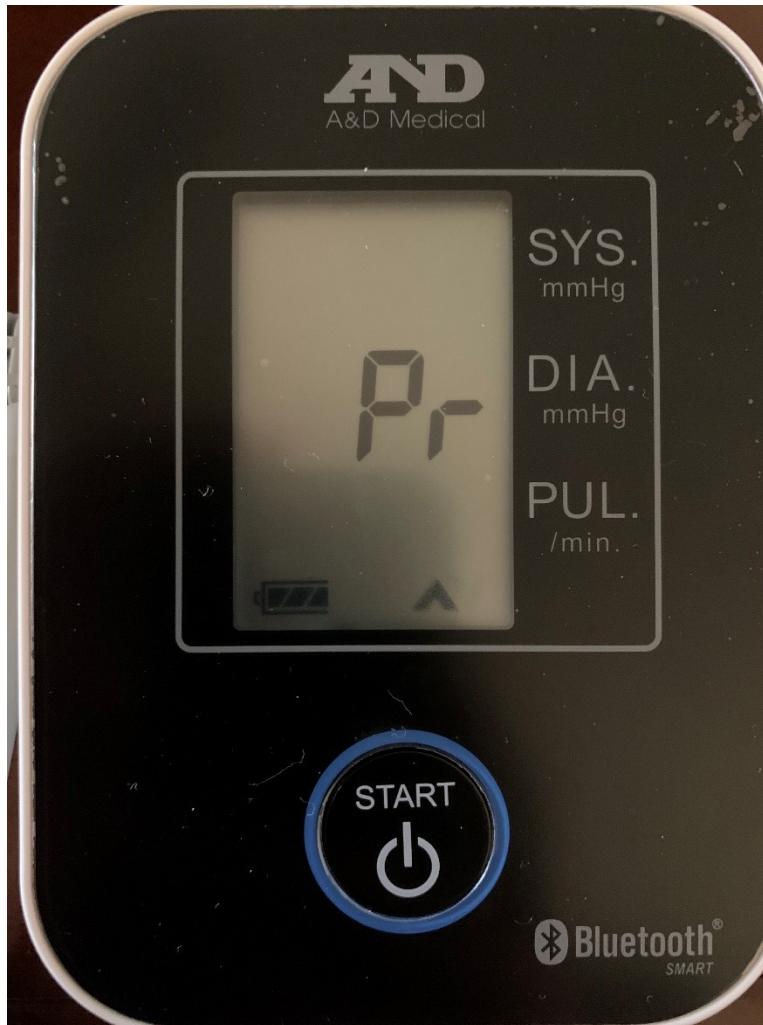

Hold down the Start button and device will go into pairing mode. Look for it to display “Pr”

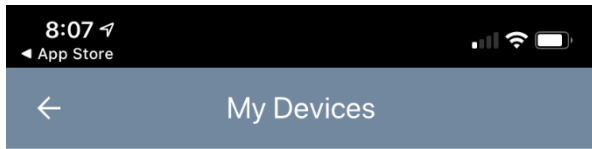


You have no devices. To add a device tap 'Add New' button.

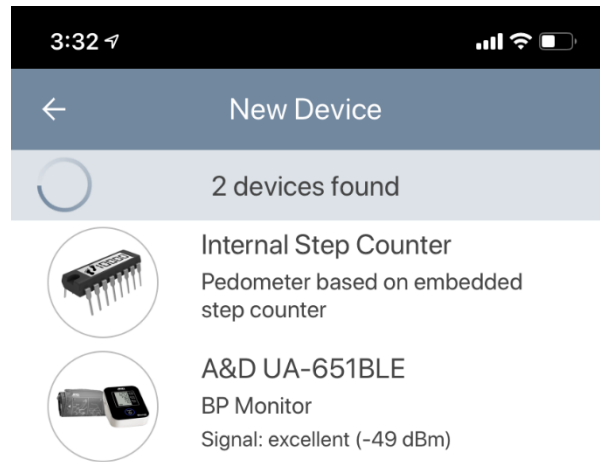


Select "Add New"



Select the A&D UA-651BLE BP Monitor

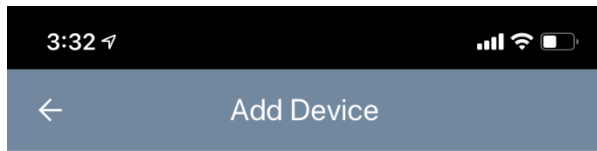

A&D UA-651BLE

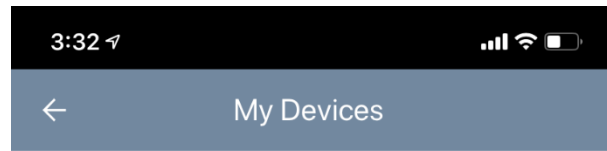

Jumper Medical
JPD-500D
Last synced at 9/28/20, 8:09 AM



A&D
UA-651BLE
Not synced yet




Select “Add to My Devices”




Back out to the Home screen




Patient App Set-Up MedM

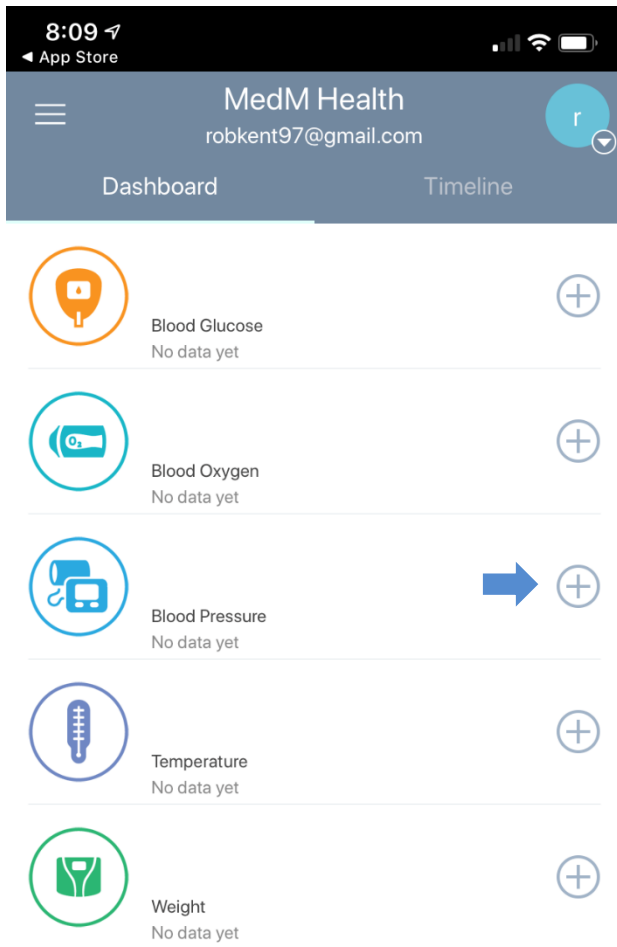

Tap the + sign next to any vital to take a reading

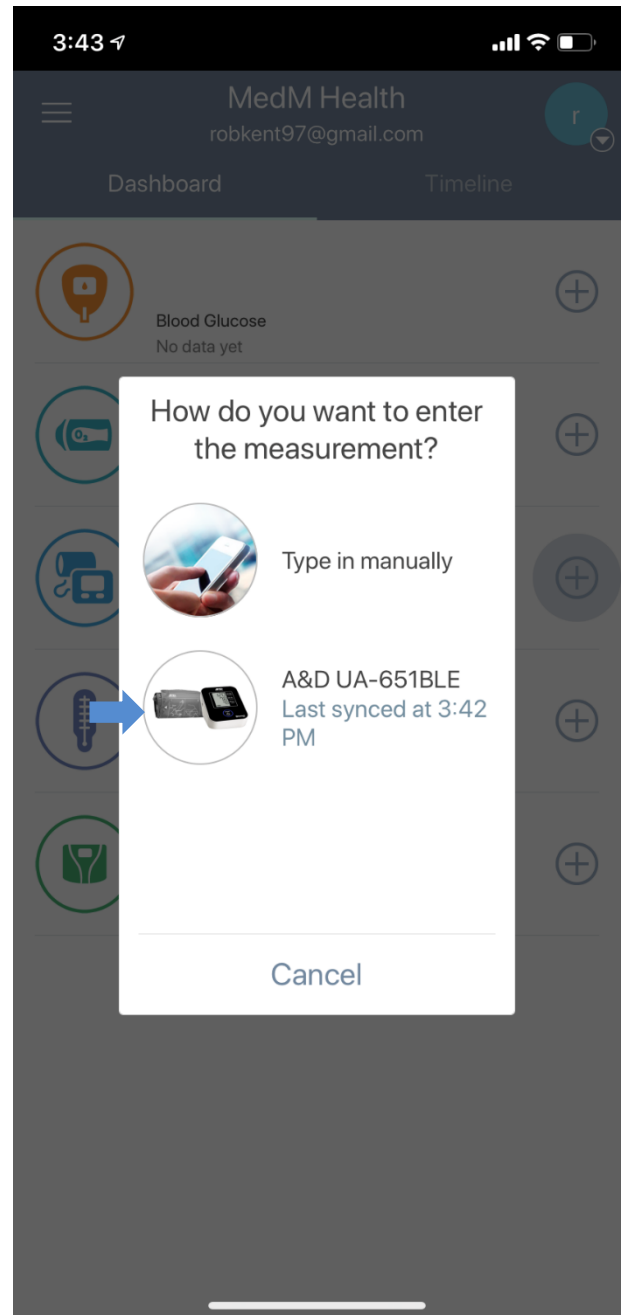

Select the synced device or manual entry

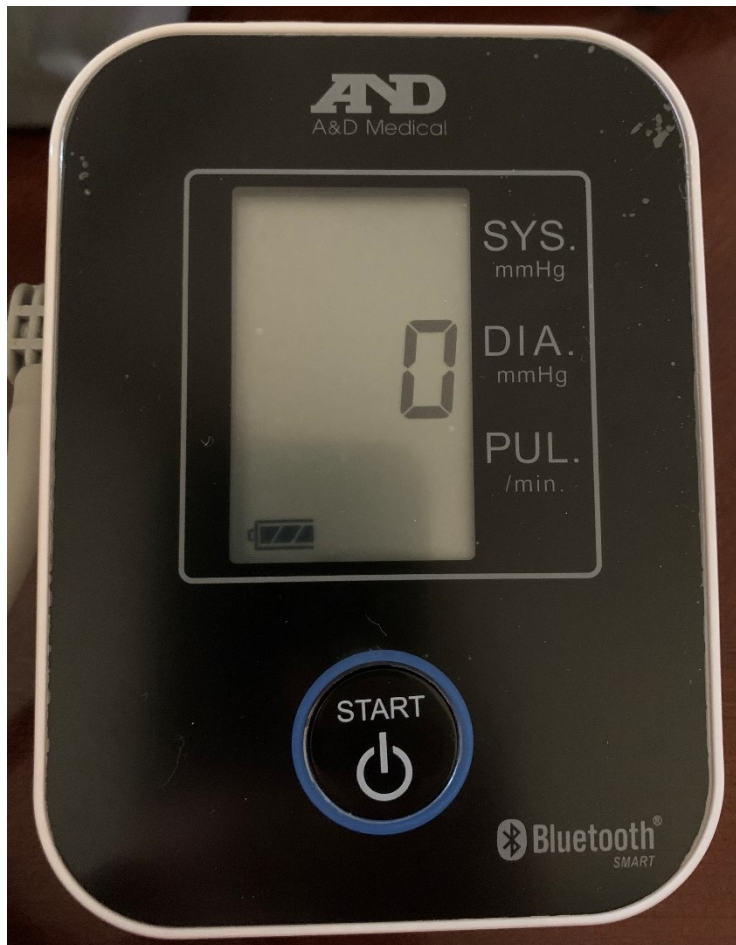

Insert arm in the cuff and secure firmly around the upper arm. Press the start button and allow the machine to cycle through a reading.



A&D UA-651BLE

Connecting to the device...



App will try to connect to the device. It will not connect until the reading is complete.



✓ 1 entry selected


✓ 27 Oct 3:44 PM ● 141/80 (43) >




Once reading is complete it will display an entry. Save the entry. Note: Device will store older readings and it may display several before or after reading is taken.




A&D Weight Scale

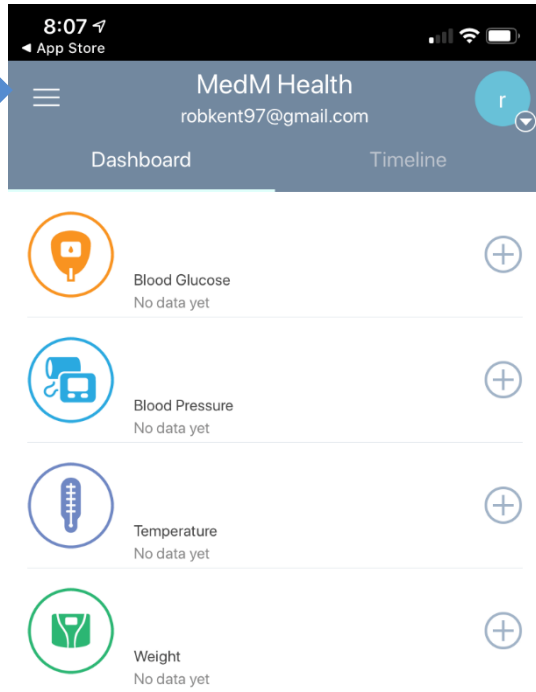


Open menu at top

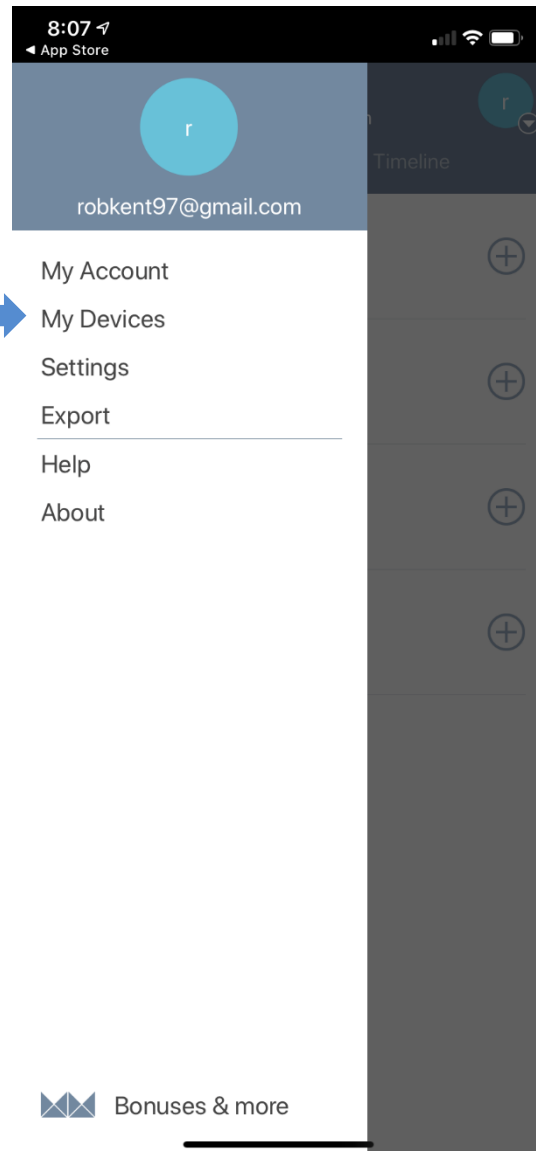


Select My Devices

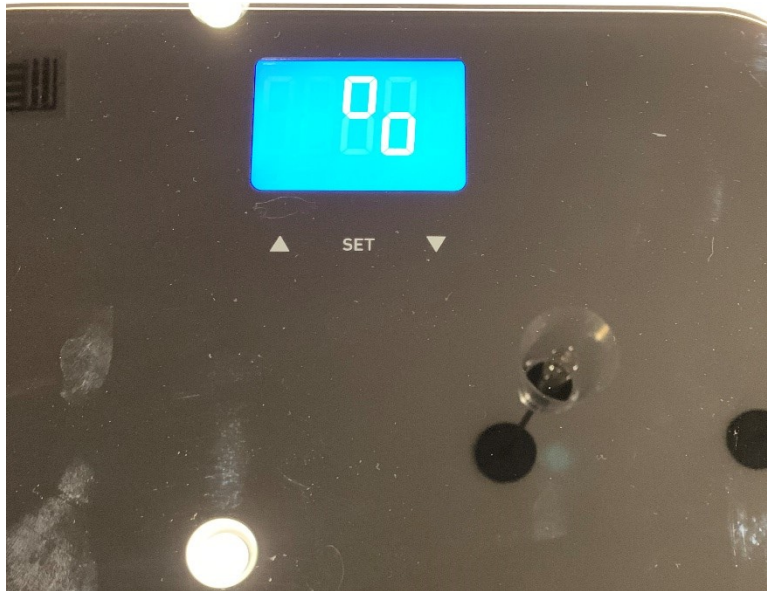

Hold down the button on the back of the scale. The display will show rotating boxes when it enters pairing mode.

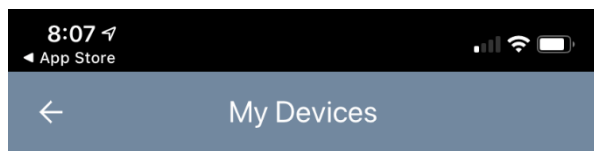




You have no devices. To add a device tap 'Add New' button.

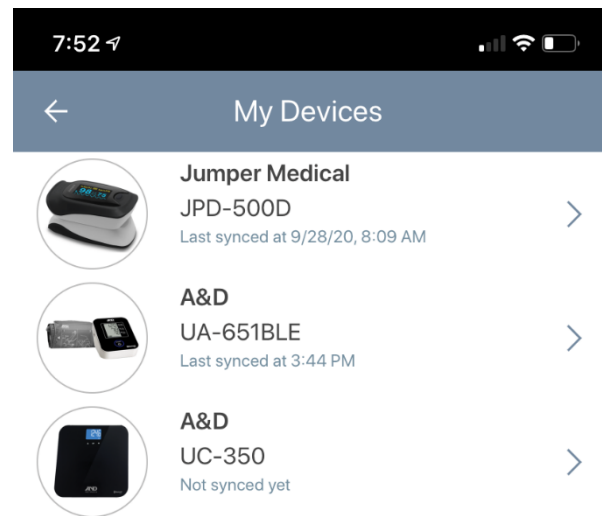

Select “Add New”




Select the A&D UC-350 Weight Scale

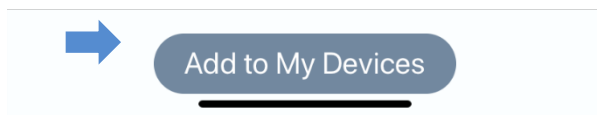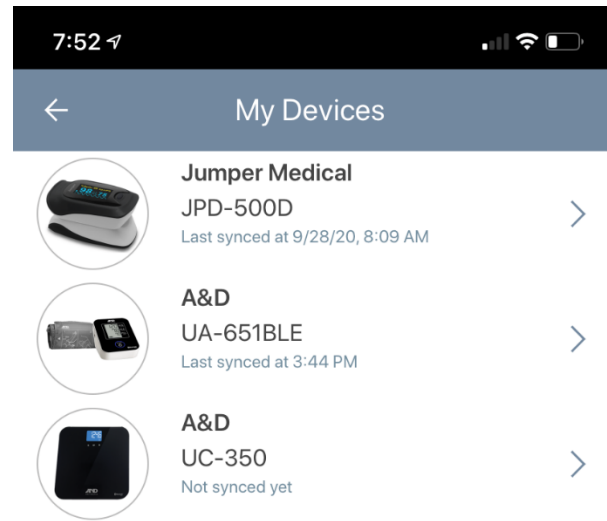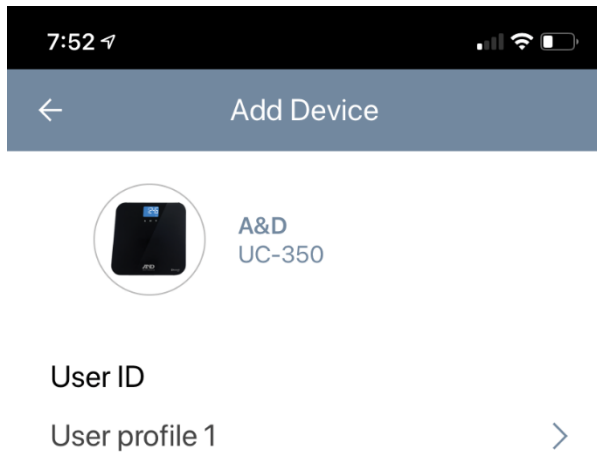

Leave on User profile 1 and Select
"Add to My Devices"




Back out to the Home screen

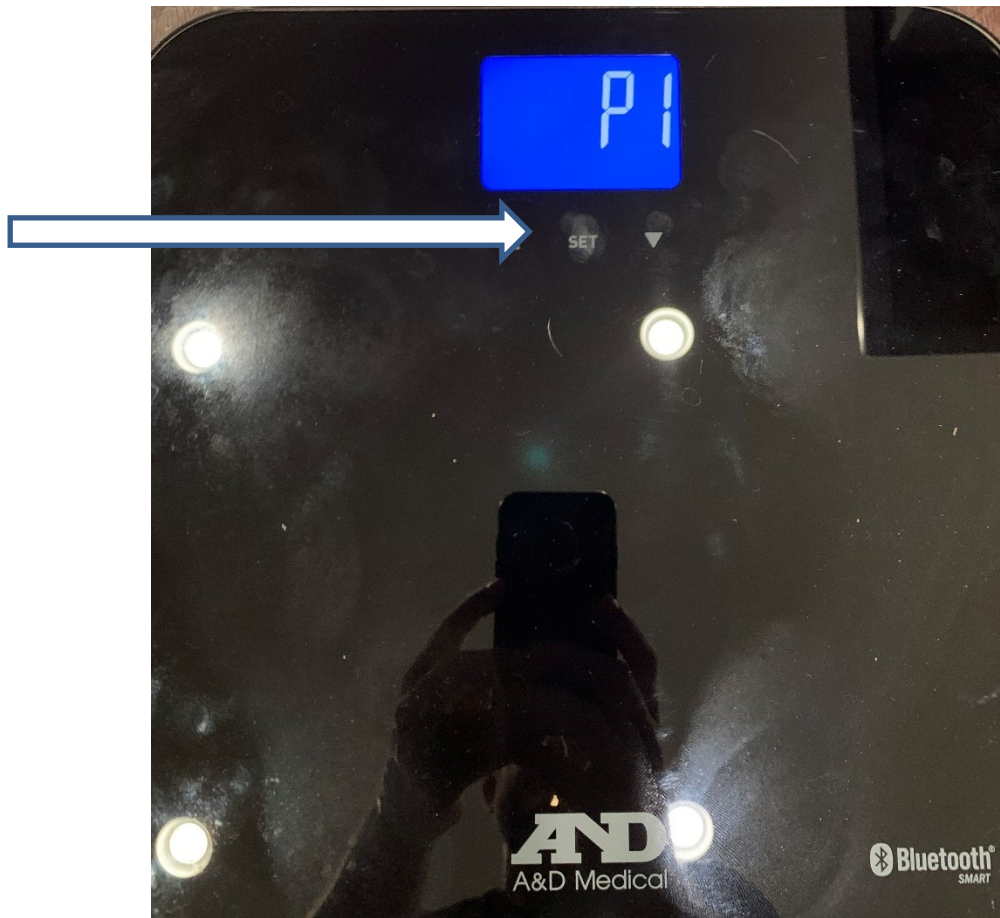


Tap the set button under the display to register Profile 1 (P1). Then all new weight readings will go to P1. Step on the scale and take a reading. The scale will save the reading and the time stamp.



Patient App Set-Up MedM

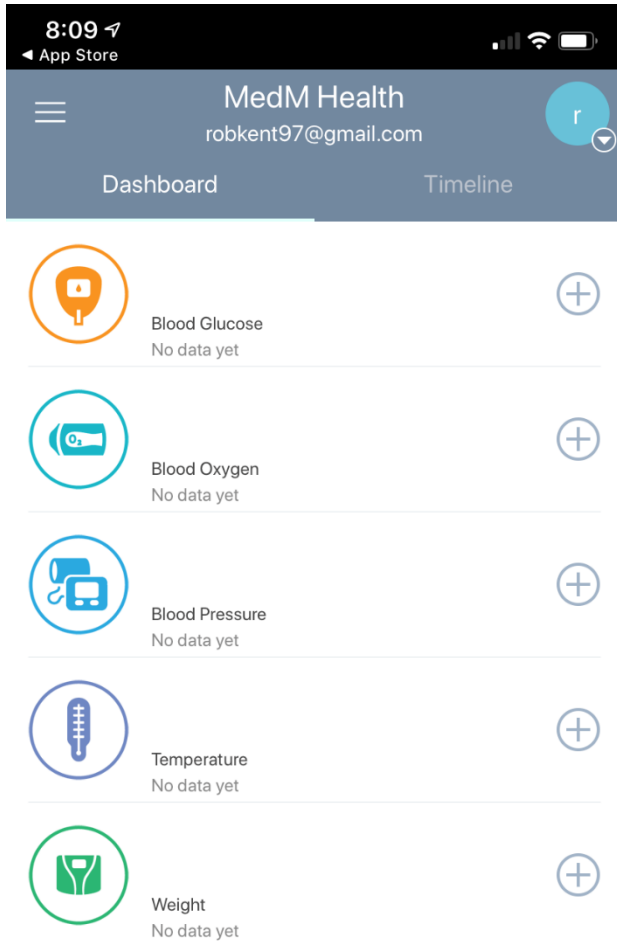

Tap the weight bar to synch with the scale and pick up any new readings. (do not use plus sign)

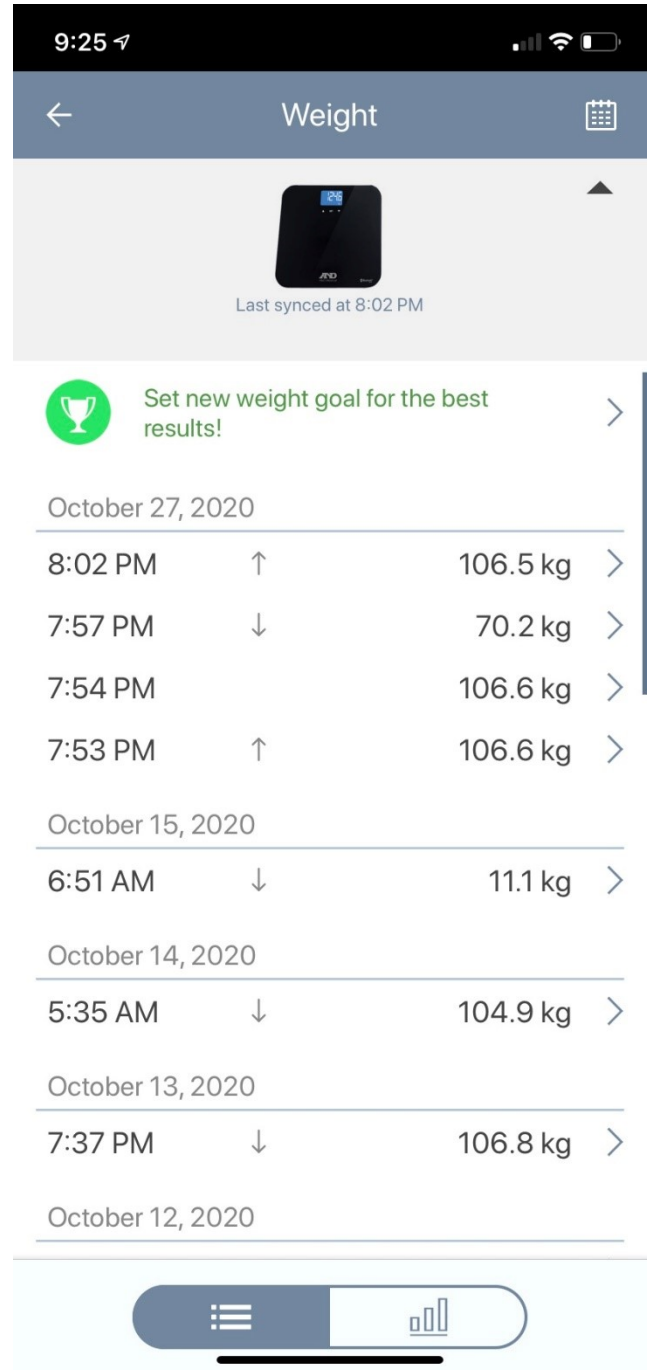

App will pick up latest reading and any other readings taken while the app was off.